


POLYMER QUANTIFICATION USING THE ROCK-EVAL® DEVICE: A METHODOLOGY FOR THE IDENTIFICATION OF PLASTICS IN SEDIMENTS



Written on 06 July 2022

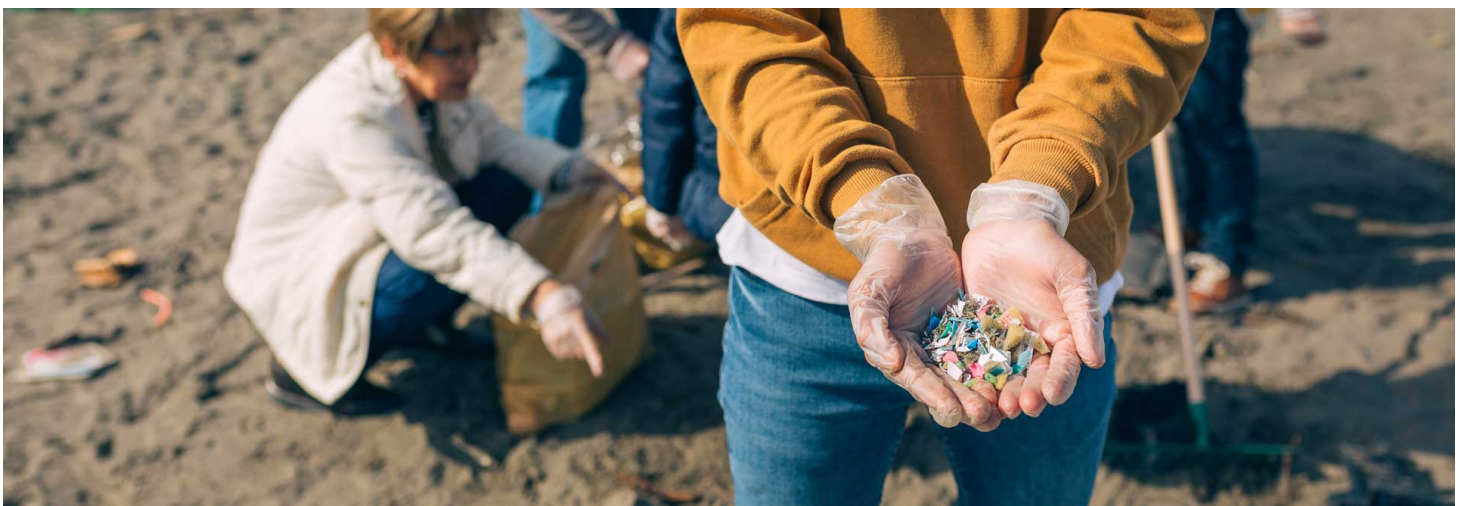

3 minutes of reading



News

Fundamental Research

Geochemistry

Chemical analysis

The presence of plastics in the ground, waterways and the oceans is of increasing concern. During the last ten years, there has been work to identify and quantify this pollution. A standardisation effort is currently underway with the objective of increasing the speed with which samples are analysed and to limit handling. IFPEN Ressources Energétiques Carnot is proposing a solution based on the Rock-Eval® device, and the initial results have just been published¹.

Rapid analysis of plastics using Rock-Eval®

The presence of plastics in the environment is a major social and economic issue^a. Significant progress has been made concerning the methods used to identify and quantify plastics in the environment over the last ten years. Nevertheless, the consensus within the scientific community that is interested in this form of pollution is that there is a need to intensify the research in order to be able to implement a standardised and rapid method that limits sample handling. IFPEN has therefore begun work on the problem by testing the potential of the Rock-Eval® device.

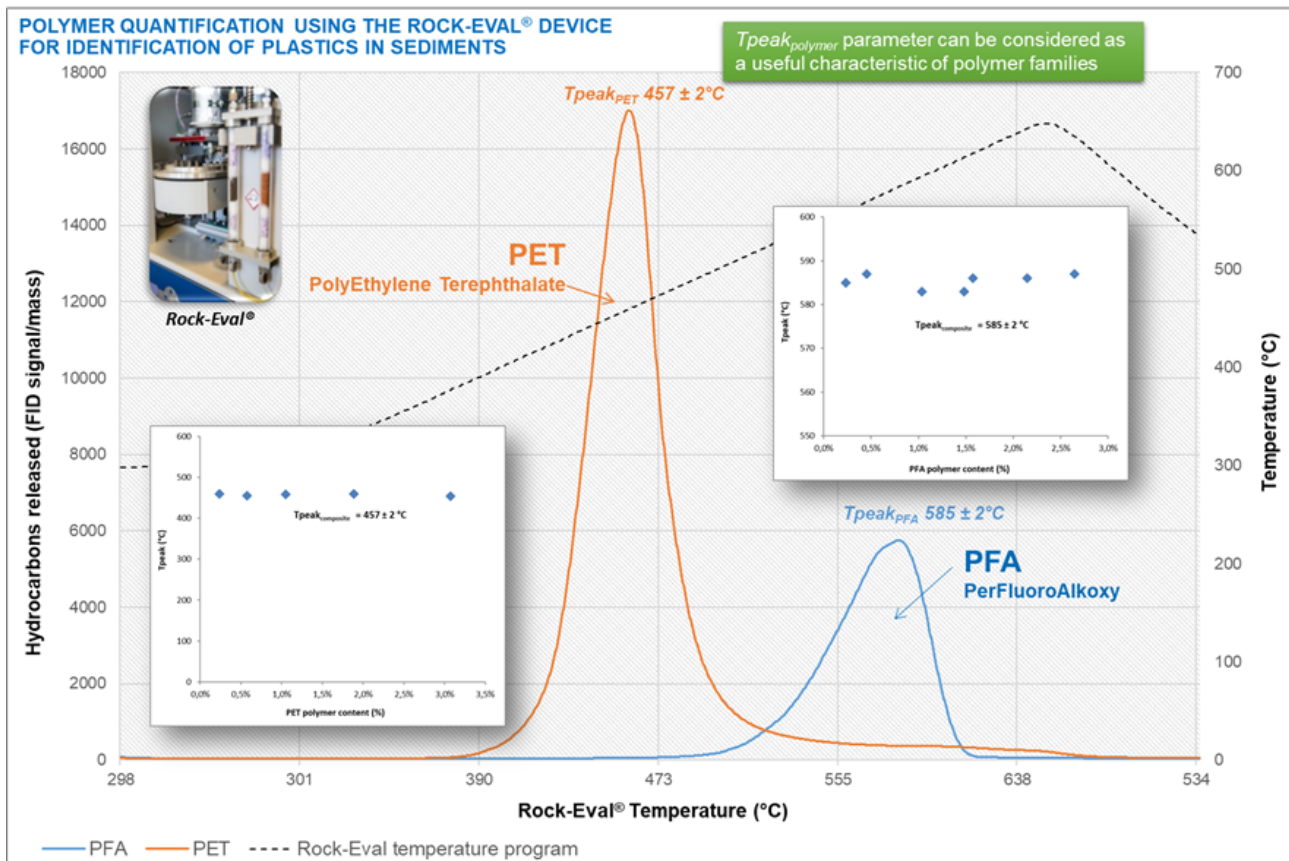


Image: Rock-Eval® thermograms illustrating the Tpeak for PET and PFA polymers

Quantification by thermal degradation of the main polymers that may contaminate the environment

This study¹ is based on synthetic samples that were previously prepared using four different mineral matrices and selected polymers (PE, PP, PE100, PA6, PA11, PFA and PET), distributed in controlled concentrations. It is demonstrated that each polymer analysed presents Rock-Eval® parameters (e.g. Tpeak^b, Total HC^c) that can be used as **specific properties** when **identifying and characterising the polymer family**. The correlation coefficients obtained for the analysed parameters and the thermal values were very good, indicating that all of these Rock-Eval® parameters can be used to define the polymer content. This implies that this thermal degradation method can be used to **characterise the plastic content in sediment samples**.

An environnement barometer capable of very precise detection

The first obtained results demonstrate **a detection limit of approximately 0.2 to 0.5 % by weight for a 60 mg sediment sample**. Additional work is being performed on this detection limit in order to extend the application of the method to as many environmental matrices as possible. The method can thus be considered **an additional and rapid screening approach** that can be implemented before or at the same time as more sophisticated analytical techniques such as SEM, Raman and FTIR measurements.

a- see for example the [2020 Senate report](#).

b- TPeak is the temperature corresponding to the maximum achieved by the flame ionization detector signal.

c- Total HC is the quantity of hydrocarbon compounds released during the total pyrolysis of the polymer.

Reference :

(1) M.-F. Romero-Sarmiento, H. Ravelojaona, D.Pillot, S. Rohais, "Polymer quantification using the Rock-Eval® device for identification of plastics in sediments", Sciences of the Total Environment, 2022. DOI : <https://doi.org/10.1016/j.scitotenv.2021.151068>

Scientific contacts : [Maria-Fermanda Romero-Sarmiento](mailto:maria-fermanda.romero-sarmiento@ifpen.fr), sebastien.rohais@ifpen.fr

YOU MAY ALSO BE INTERESTED IN

[Rock-Eval® : Thermal analysis of rocks and soils](#)

[Rock-Eval®, supporting soil research for the climate challenge](#)

[Rock-Eval® : to find out more](#)

Link to the web page : [Polymer quantification using the Rock-Eval® device: a methodology for the identification of plastics in sediments](#)